


Radio Frequency (RF) microneedling system use specially insulated needles that deliver high-intensity radio frequency energy into the targeted tissue once the needles have broken through the skin. This gives a uniform warming effect to the deep tissue, which helps to stimulate new collagen fibers to grow, leading to significant improvement to skin quality and texture.

Therapy system RF Fractional Microneedle machine for wrinkle removal and skin tightening.

Technology

Israel's most advanced radio frequency technology, use matrix needles to penetrate 0. 5mm-3.5mm subcutaneously, and at the same time, each pulse emits radio frequency energy through multiple matrix points, which is directly heated under the skin; Under the stimulation of matrix needles, subcutaneous produce a large amount of collagen, and use radio frequency to directly heat under the skin to repair broken collagen fibers, so as to achieve a better repair effect.

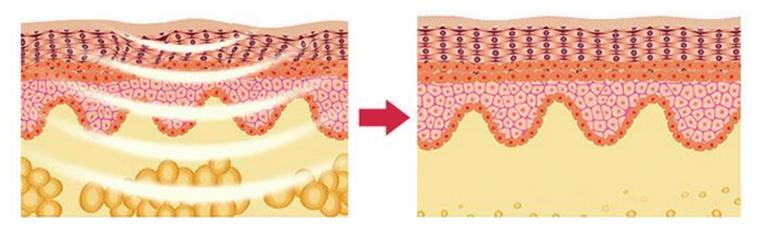

Israel's most advanced radio frequency technology

Principle

Microneedling with radiofrequency, like other microneedling procedures, creates controlled micro-injuries in the skin to stimulate collagen production. Other methods of microneedling utilize PRP, stem cell growth factor, and other serums to penetrate deep into the skin and further enhance results.

After cleansing the skin surface and applying a topical anesthetic, the microneedling device is gently pressed against the treatment area to create numerous, tiny, microchannels. Microneedling with radiofrequency uses the added technology of flow needles to deliver RF energy into the dermis. The RF energy heats the dermis which not only encourages collagen growth but tissue tightening as well.

The penetration of the microneedles into the skin leads to a release of growth factors and triggers the wound healing cascade, encouraging collagen production which leads to younger looking skin. The needles also help to mechanically break down the scar tissue. Since the epidermis is not destroyed, as is the case with all microneedling treatments, recovery time is much shorter when compared to more aggressive laser resurfacing treatments or deep chemical peels.

The goal of a microneedling with radiofrequency procedure is to target imperfections of the skin with micro-injuries which cause the release of growth factors and triggers the wound healing process which produces collagen, creating younger, tighter looking skin. The stamping action of the Microneedling RF Therapy System's vertical tip creates invisible micro-injuries to the skin, while stimulating the skin's natural ability to heal itself and increase the production of new collagen. In doing so, these micro-injuries encourage and harness the body's innate ability to regrow and repair the skin through normal physiological processes and the skin becomes rejuvenated and refreshed.

LA REINA



Activating the body's natural healing processes by applying RF assisted ultra-thin electrodes ablation, followed by galvanic current to the tissue to create a thermal effect creating micro channels to induce controlled micro wounds, thus encouraging the generation of new and healthy skin.

Selectable tip Type

The PrettyLasers's RF micro-needling device provides a unique combination of treatment heads that let you treat a wide range of applications for all skin types, including challenging skin.

1 suitable for eyes

2 suitable for forehead

3 & 4 suitable for face and neck

PLUMA

Applying contact RF and pin point ablation with controlled volumetric heating to achieve skin resurfacing and regeneration via wound healing process, resulting in skin texture improvement and refinement.



Adjustable Needle Depth


PrettyLasers's RF micro-needling features adjustable needle penetration ranging from 1.5-3.0 mm offering customized dermal remodeling based on the area and the depth of the condition being treated.


PrettyLaser's Microneedling RF System treats to the maximum depth of 3.0mm – which is the deepest amongst RF microneedling technology. This means that even the deepest acne scars can be treated effectively.



User-friendly Interface

Graphical User Interface - It is easy to operate and learn.

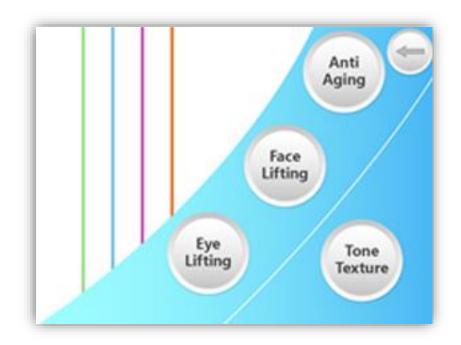

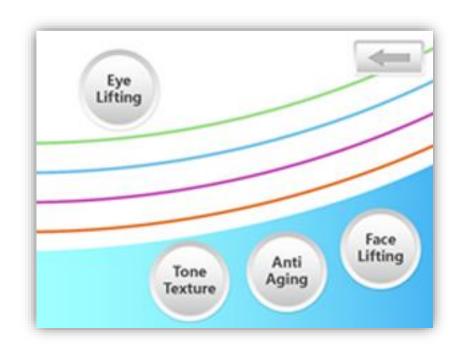

Technical Parameters

Power	300W
Treatment Area	13mm (round head), 17mm*33mm
Frequency	2Mhz
Screen	15inch Touch Screen
Voltage	AC220V, 50/60HZ
Working Environment	Ambient temperature 25-35 degree, humidity <95%
Package size	52*65*118cm
Head	10 pcs changeable heads
Smooth heads	4 pcs
25 pins Heads with needles	25 pins with size 1.5mm 2.0mm 2.5mm
49 pins Heads with needles	49 pins with size 2.0mm 2.5mm 3.0mm

Professional Microneedling Therapy System Benefits

The gold standard for skin rejuvenation

Whitening and rejuvenating skin tighten pores, brighten skin.

→ Diminish wrinkles & lines

Smooth stretch marks

→ Diminish acne scars and other scarring

Whatsapp: 0086 15001300057

Web: https://www.prettylasers.com/

Email: info@prettylasers.com

Phone: +0086 15001300057



Contact With PrettyLasers